

APPLICATION FOR OSHPD SPECIAL SEISMIC	OFFIC	E USE ONLY
CERTIFICATION PREAPPROVAL (OSP)	APPLICATION #:	OSP – 0564 – 10
OSHPD Special Seismic Certification Preapproval (OSP)		
Type: 🛛 New 🗌 Renewal		
Manufacturer Information		
Manufacturer: Siemens Healthcare GmbH, Diagnostic Imaging, Com	puted Tomography	
Manufacturer's Technical Representative: Ottmar Förstel		_
Mailing Address: Siemensstr. 3, 91301 Forchheim, Germany		
Telephone:	.foerstel@siemens.com	<u>1</u>
Product Information FOR CODE	01	
Product Name: SOMATOM CT Systems	The states	
Product Type: Computed Tomography (CT) medical imaging system	TE.	
Product Model Number: See Attachment	G	
(List all unique product identification numbers and/or part numbers)		
General Description: <u>Multiple component system for producing Con</u>	nputed Tomography (C	T) medical images for a
wide variety of medical diagnostic results.	0	
Mounting Description: Rigid floor mounted, see attachment.		
	(V	
Applicant Information	CODÉ	
Applicant Company Name:		
Contact Person: Travis Soppe, SE		
Mailing Address: _ 250 Bobwhite Ct, Suite 100, Boise, ID 83706		
Telephone: (208) 342-5898 Ext. 115 Email: tsoppe	@wegai.com	
I hereby agree to reimburse the Office of Statewide Health I accordance with the California Administrative Code, 2016.	Planning and Devel	opment review fees in
Signature of Applicant:	Date	e: 05-07-2018
Title: President Company Name: W.E. C	Gundy & Associates, Inc	2
"Access to Safe, Quality Healthcare Environments that Meet California's Diverse and Dynamic Needs"	A LA AM JAAA	OSHPD

"Access to Safe, Quality Healthcare Environments that Meet California's Diverse and Dynamic Needs" STATE OF CALIFORNIA – HEALTH AND HUMAN SERVICES AGENCY OSH-FD-759 (REV 12/16/15)



California Licensed Structural Engineer Responsible for the Engineering and Test Report(s)
Company Name:W.E. Gundy & Associates, Inc.
Name: Travis Soppe, SE California License Number: S6115
Mailing Address: _250 Bobwhite Ct, Suite 100, Boise, ID 83706
Telephone: (208) 342-5898 Ext. 115 Email: tsoppe@wegai.com
Supports and Attachments Preapproval
Supports and attachments are preapproved under OPM- (Separate application for OSHPD Preapproval of Manufacturer's Certification (OPM) of Supports and attachments is required)
Supports and attachments are not preapproved CODR
Certification Method
Testing in accordance with: ICC-ES AC156 Other (Please Specify): OSP-0564-10
Testing Laboratory
DATE: 02/21/2019
Company Name: IABG mbH Contact Name: Dr. Steffen Roedling
Mailing Address: Einsteinstrasse 20, Ottobrunn, Germany D-85521
Telephone:

"Access to Safe, Quality H	lealthcare Environments th	at Meet California's	Diverse and	Dynamic Need	s
----------------------------	----------------------------	----------------------	-------------	--------------	---

OFFICE OF STATEWIDE HEALTH PLANNING AND DEVELOPMENT FACILITIES DEVELOPMENT DIVISION

Seismic	Parameters
---------	------------

Design in accordance with ASCE 7-10 Chapter 13: 🛛 Yes 🗌 No
Design Basis of Equipment or Components (Fp/Wp) = <u>See Attachment</u>
S_{DS} (Design spectral response acceleration at short period, g) = 2.00 (z/h = 1); 2.50 (z/h = 0)
a _p (In-structure equipment or component amplification factor) = <u>See attachment</u>
R _P (Equipment or component response modification factor) = <u>See attachment</u>
Ω_0 (System overstrength factor) = _See attachment
I _p (Importance factor) = 1.5
z/h (Height factor ratio) = _1.0 (S _{DS} = 2.00); 0 (S _{DS} = 2.50)
Equipment or Component Natural Frequencies (Hz) = See attachment
Overall dimensions and weight (or range thereof) = See attachment
Equipment or Components @ grade designed in accordance with ASCE 7-10 Chapter 15: Yes No Design Basis of Equipment or Components (V/W) =
S _{DS} (Design spectral response acceleration at short period, g) =
S_{D1} (Design spectral response acceleration at 1 second period, g) = $\sqrt{2}$
R (Response modification coefficient) = <u>OSP-0564-10</u>
Ω_0 (System overstrength factor) =
C₄ (Deflection amplification factor) = ^{BY} :AL1 Sumer
I_p (Importance factor) = 1.5
Height to Center of Gravity above base =
Equipment or Component Natural Frequencies (Hz) =
Overall dimensions and weight (or range thereof) =
Tank(s) designed in accordance with ASME BPVC, 2015: 🔲 Yes 🖾 No
List of Attachments Supporting Special Seismic Certification
🛛 Test Report(s) 🗌 Drawings 📄 Calculations 🖾 Manufacturer's Catalog
 ✓ Other(s) (Please Specify): Certified System Matrix, UUT Summary Sheets, Subcomponent Certification Letter
OSHPD Approval (For Office Use Only) – Approval Expires on December 31, 2022
Signature: Date: February 20, 2019
Print Name: Ali Sumer Title: DSE
Special Seismic Certification Valid Up to : S _{DS} (g) = <u>See Above</u> z/h = <u>See Above</u>
Condition of Approval (if applicable):
"Access to Safe, Quality Healthcare Environments that Meet California's Diverse and Dynamic Needs"
STATE OF CALIFORNIA – HEALTH AND HUMAN SERVICES AGENCY OSH-FD-759 (REV 12/16/15) Page 3 of 3

SIEMENS HEALTHCARE GmbH SPECIAL SEISMIC CERTIFICATION CERTIFIED SYSTEM AND COMPONENTS



Manufacturer: Siemens Healthcare GmbH

System: SOMATOM CT Systems

	Siemens]	Dimensions (in)			TITIT
System Component ¹	Part Number	Width	Length	Height	- Weight (lb)	Mounting	UUT
		FOBGant	ries				
SOMATOM Edge Plus	10267000	93.7	D 36.8	78.0	4607	floor	UUT-1
SOMATOM Definition Edge	8098555	93.7	36.8	78.0	4851	floor	interpolated
SOMATOM Force	10742326	94.5	46.5	F178.3	5900	floor	UUT-2
·		age Reconstr	uction Syster	ns			
PC IRSmx5a	10590110	8.3 DATE: 02/2	17.7	<mark>_2</mark> 6.8	52.8	floor	UUT-5
PC IRSmx5c	10590112	8.3	17.7	26.8	60.0	floor	interpolated
PC IRSmx5d	10590113	8.3	17.7	26.8	60.0	floor	interpolated
PC IRSmx5e	10590114	A 8.3	17.7	26.8	60.0	floor	interpolated
PC IRSmx5b	10590111	8.3	17.7	26.8	61.6	floor	UUT-6

Notes:

¹ All components are manufactured by Siemens Healthcare GmbH unless noted. Part numbers listed uniquely identify type of component, manufacturer, and material of construction for each sub-componenent within the tested units.

SIEMENS HEALTHCARE GmbH SPECIAL SEISMIC CERTIFICATION CERTIFIED SYSTEM AND COMPONENTS



Manufacturer: Siemens Healthcare GmbH

System: SOMATOM CT Systems

	Siemens	Dimensions (in)				Weight (
System Component ¹	Part Number	Width	Length	Height	Table	Patient	Total	Mounting	UUT
			PHS Patie	nt Tables		-			
PHS-4	8097144	29.5	95.7-158.7	18.9-36.2	933	0 - 415	933 - 1348	floor	UUT-4x
PHS-5	10742323	29.5	100.5-179.5	24.3-41.7	1066	0 - 415	1066 - 1475	floor	interpolated
PHS-4n	11268204	29.5	100.5-179.5	6 24.3 41.7	1060	0 - 415	1060 - 1475	floor	UUT-3
	· · ·	R	MPT Patie	ent Tables	E	•			
MPT-2	8097102	2 <mark>9.5</mark>	95.7-174.4	21.6-36.2	1 <mark>311</mark>	0 - 560	1311 - 1870	floor	UUT-5x
MPT-2n	11268202	29.5	D9613-17520 2	1/22.0-36.2	1175	0 - 560	1175 - 1735	floor	interpolated
MPT-4	10742324	27.8	99.6-179.5	23.0-38.6	1258	0 - 560	1258 - 1818	floor	UUT-4

Notes:

- Patient tables utilize the same active components for each PHS-x and MPT-x configurations with the only difference being the non-active cover configuration.

- The "x" at the end of the UUT designation indicates a test was performed at a different time than the other UUT's.

- The UUT summary sheets detail the normal operating position used for the test and the simulated patient weight.

¹All components are manufactured by Siemens Healthcare GmbH unless noted. Part numbers listed uniquely identify type of component, manufacturer, and material of construction for each sub-componenent within the tested units.

TABLE 3

SIEMENS HEALTHCARE GmbH SPECIAL SEISMIC CERTIFICATION CERTIFIED SUBCOMPONENTS



Manufacturer: Siemens Healthcare GmbH

System: SOMATOM CT Systems

C L 	Siemens	D	imensions (Weight	UUT	
Subcomponent ¹	Part Number	Width	Length	Height	(lb)	UUI
		Collimat	or			
Type Z80S A	10742700	18.1	10.7	5.5	42	UUT-2
Type Z80S B	10742701	18.1	10.7	5.5	42	UUT-2
Type Z65	10589900	18.1	10.4	5.1	36	UUT-1
	X-Ra	y Tube A	ssembly			
2x Vectron	10414460	R 16.0D	16.9	11.1	144	UUT-2
Straton MX Sigma	11270277	17.0	9.8	11.1	84	UUT-1
Straton MX	8401825	17.0	9.8	11.1	84	interpolated
	Hig	h Voltag	e Tank	T.		
HVT MxT (2x)	10513610 C	SP- 1958 64	-1012.3	12.4	89	UUT-2
HVT Matrix	<mark>836</mark> 5707	11.8	19.3	12.4	89	UUT-1
·	O BY:A	Detecto	per	0		
DMS P58A	10393130	45.7	14.1	17.8	154	UUT-2
DMS P58B	10393140	32.4	13.3	15.8	107	UUT-2
DMS P46F	10186150	46.0	18.7	12.2	151	UUT-1
	PV SI	ipring Sy	stem	ès /		
Type P58	10494903	ALL 3.0	62.2	62.2	152	UUT-2
Type P46	10187169	2.8	59.8	59.8	187	UUT-1
	С	ooling Sy	stem			
Water Cooling - Glen Dimplex -	10742370	20.0	30.0	21.7	143	UUT-2
Water Cooling - Glen Dimplex -	10589646	16.7	25.7	20.3	130	UUT-1

¹ All components are manufactured by Siemens Healthcare GmbH unless noted below the subcomponent name.

UNIT UNDER TEST (UUT) SUMMARY SHEET




Mounting Details: Rigid floor mounted with 4 - M16 grade 12.9 bolts.



UNIT UNDER TEST (UUT) SUMMARY SHEET




Mounting Details: Rigid floor mounted with 4 - M16 grade 12.9 bolts.

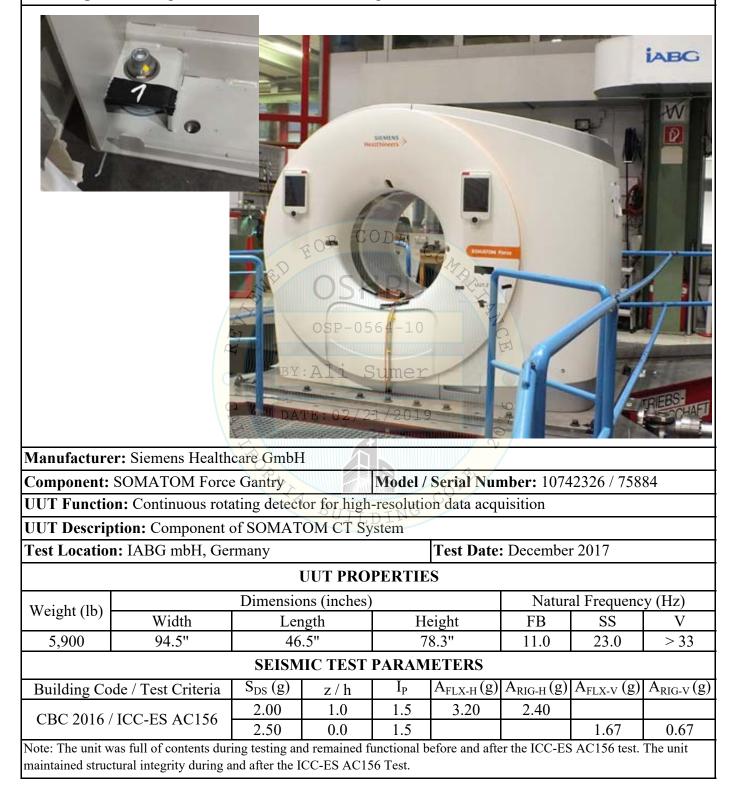

UNIT UNDER TEST (UUT) SUMMARY SHEET



Mounting Details: Rigid Floor mounted with 4 - M10 grade 12.9 bolts

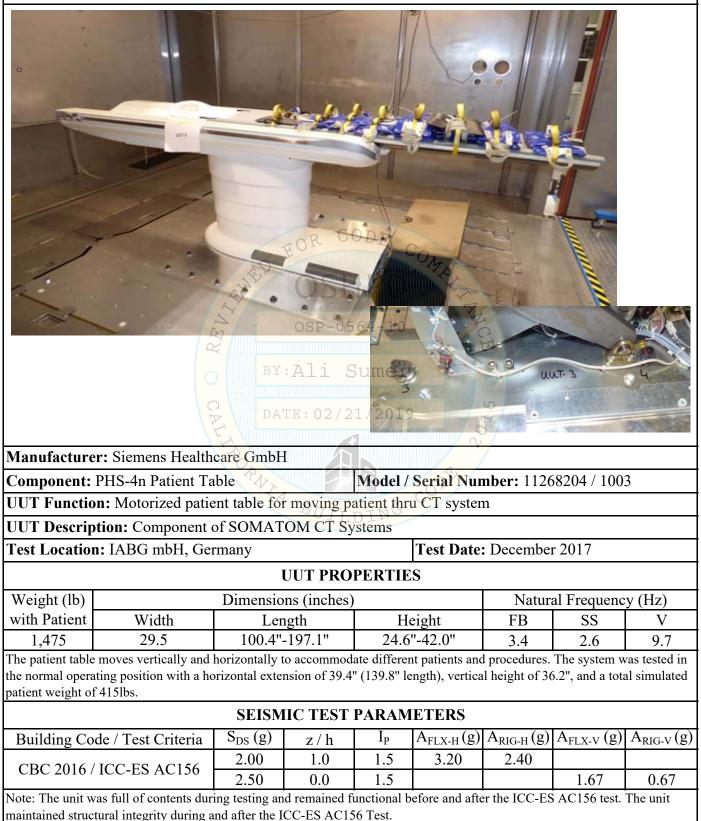

UNIT UNDER TEST (UUT) SUMMARY SHEET



Mounting Details: Rigid Floor mounted with 4 - M10 grade 12.9 bolts



UUT-4x

UNIT UNDER TEST (UUT) SUMMARY SHEET



Mounting Details: Rigid floor mounted with 4 - M10 bolts



UUT-5x

UNIT UNDER TEST (UUT) SUMMARY SHEET



Mounting Details: Rigid floor mounted with 4 - M10 bolts




UUT-5

UNIT UNDER TEST (UUT) SUMMARY SHEET



Mounting Details: Rigid Floor mounting using Siemens provided seismic restraint kit SN:10432402. Seismic restraint kit includes a 1" wide hand tightened cam buckle strap (560lb WLL) looped thru angle brackets positioned on each side of the unit. The angle brackets are attached to the table with individual M10 grade 8.8 bolts.

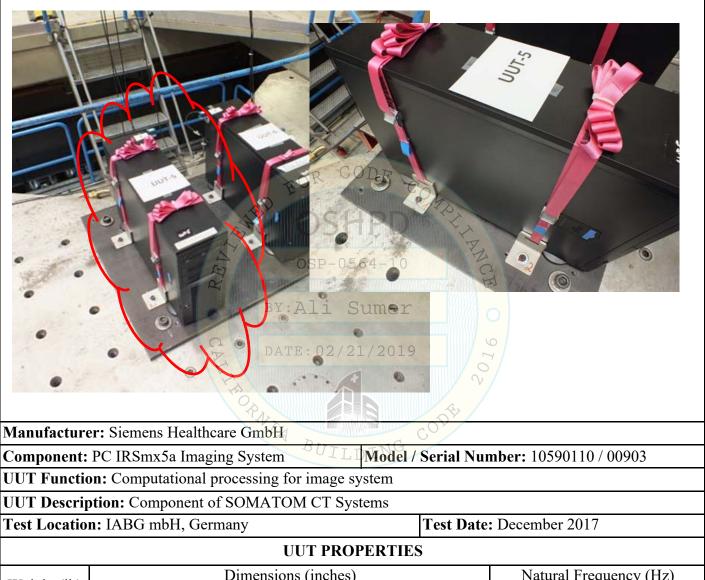

Weight (lb)		Dimensions (inches)						Natural Frequency (Hz)			
weight (10)	Width	De	Depth		Height		SS	V			
53	8.3"	17	.7"	2	6.8"	> 33	27.4	> 33			
	SEISMIC TEST PARAMETERS										
Building C	ode / Test Criteria	$S_{DS}(g)$	z / h	I _P	$A_{FLX-H}(g)$	$A_{RIG-H}(g)$	$A_{FLX-V}(g)$	$A_{RIG-V}(g)$			
CPC 2016	/ ICC-ES AC156	2.00	1.0	1.5	3.20	2.40					
CBC 2010	/ ICC-ES AC150	2.50	0.0	1.5			1.67	0.67			
Note: The unit w	vas full of contents during	g testing and	remained fur	ctional bef	ore and after	the ICC-ES A	AC156 test. T	he unit			


Note: The unit was full of contents during testing and remained functional before and after the ICC-ES AC156 test. The unit maintained structural integrity during and after the ICC-ES AC156 Test.

UUT-6

UNIT UNDER TEST (UUT) SUMMARY SHEET



Mounting Details: Rigid Floor mounting using Siemens provided seismic restraint kit SN:10432402. Seismic restraint kit includes a 1" wide hand tightened cam buckle strap (560lb WLL) looped thru angle brackets positioned on each side of the unit. The angle brackets are attached to the table with individual M10 grade 8.8 bolts.

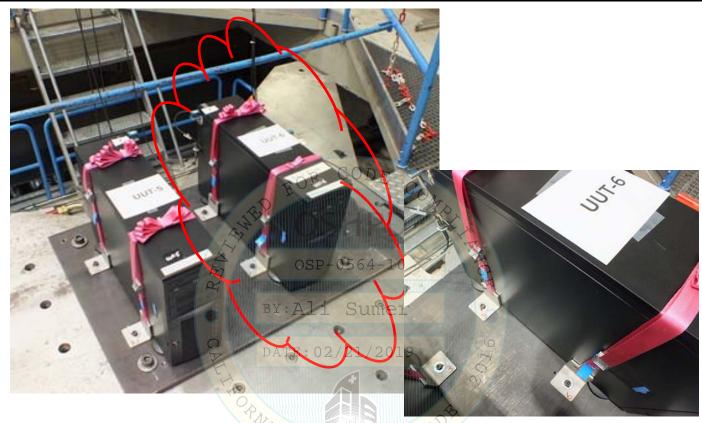

		17.						
Manufacture	er: Siemens Healthca	are GmbH	Da	(IIII)	COR			
Component:	PC IRSmx5b Imagi	ng System	POILI	Model /	Serial Nur	nber: 1059	0111 / 009	03
UUT Functio	on: Computational p	rocessing fo	or image sy	vstem				
UUT Descrip	otion: Component of	SOMATO	M CT Syst	tems				
Test Location	n: IABG mbH, Gern	nany			Test Date	: Decembe	r 2017	
		U	UT PROP	PERTIE	S			
Waight (1b)	Dimensions (i				Natural Frequency (H			y (Hz)
Weight (lb)	Width	De	pth	Н	leight	FB	SS	V
62	8.3"	17	.7"	2	26.8"	> 33	30.3	> 33
		SEISM	IC TEST I	PARAM	ETERS			
Building Co	ode / Test Criteria	$S_{DS}(g)$	z / h	I _P	$A_{FLX-H}(g)$	$A_{RIG-H}(g)$	$A_{FLX-V}(g)$	$A_{RIG-V}(g)$
CPC 2016	/ ICC-ES AC156	2.00	1.0	1.5	3.20	2.40		
CBC 2010	/ ICC-ES AC130	2.50	0.0	1.5			1.67	0.67
	vas full of contents durin tural integrity during an	0 0			fore and after	the ICC-ES A	AC156 test. T	ne unit